

# ATTP MOBILE USER GUIDE


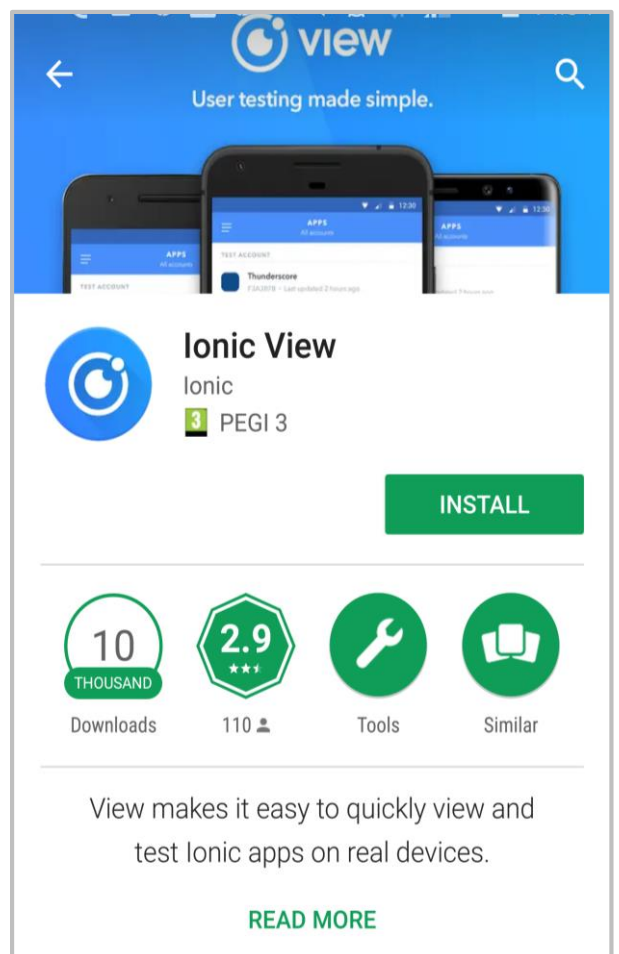
## INSTALLATION

- Install **Ionic View** from **App Store** or **Google Play**  
(not Legacy Ionic View, it is old version)



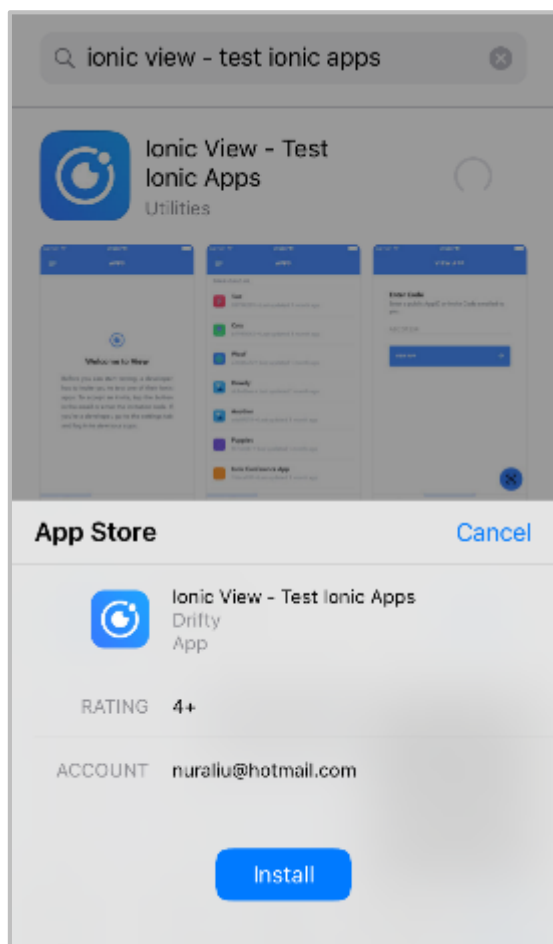
- Open **View App** tab
- Enter AppID: **957c3074**
- Click **View App**
- **ATTPMobile** app will start
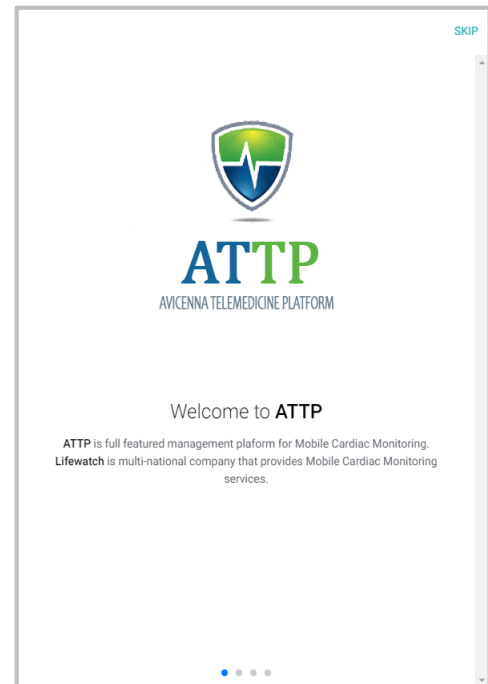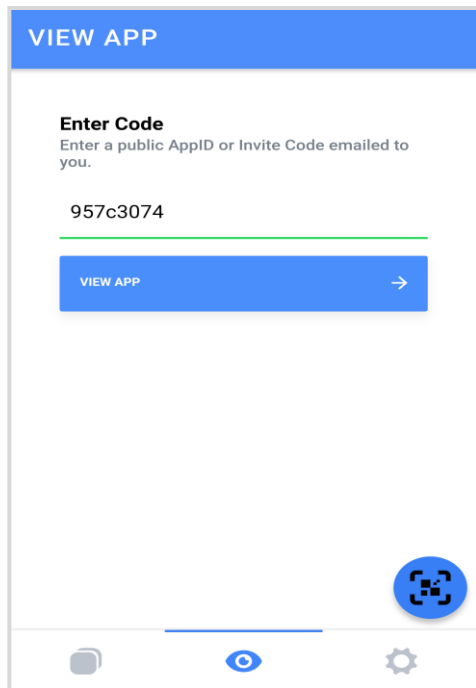
## APPLICATION USAGE

- Open **Ionic View**
- Click **ATTPMobile** from Applications list
- **ATTPMobile** will start
- Skip introduction page
- Enter your **Username & Password**  
(only accounts having **Sistem Yöneticisi** role are allowed for now)
- Select the environment (PROD)
- Click **Login**
- If authorized, you will see the list of enrollments
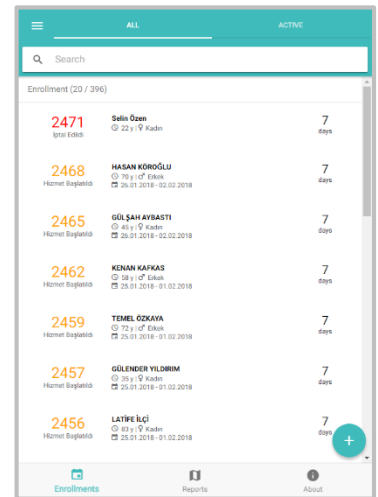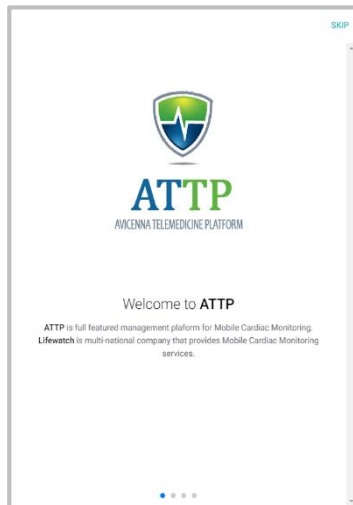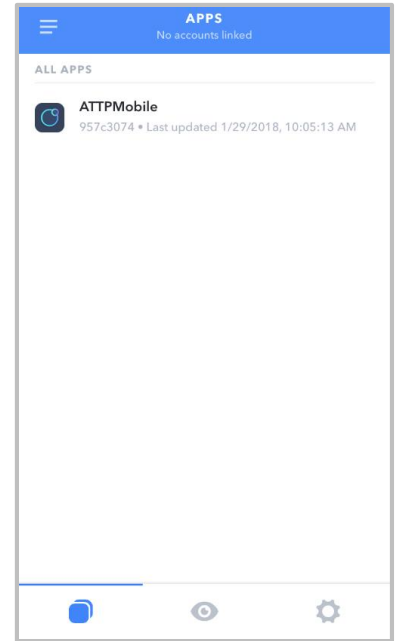
## MAJOR SCREENS

- Intro Page
- Login Page
- Enrollments List
- Reports List
- Profile Page
- About Page

## MAJOR FUNCTIONALITIES

- Multi Deployment Access
- List/Filter Enrollments
- List/Filter Reports
- Report Download
- Call Patient
- Call Doctor
- Filter by Status (Active/All)
- Filter by Keyword
- Filter by Report Type

## MAJOR FEATURES

- Material UI Design
- Infinite Scrolling (on-demand)
- Conditional Formatting
- Side Menu Access
- Deployed on both App Store & Google Play

# SCREENSHOTS

SKIP

**ATTP**  
AVICENNA TELEMEDICINE PLATFORM

## Welcome to **ATTP**

**ATTP** is full featured management platform for Mobile Cardiac Monitoring.  
**Lifewatch** is multi-national company that provides Mobile Cardiac Monitoring services.

SKIP

## Start using ATTP

You can start using ATTP just by reaching Lifewatch in Turkey.  
Then you will have an account

SKIP

## Start getting ECG Reports

After getting your **doctor account** you will be able to create new **enrollments** for your patient

And then you will start getting live ECG reports from you patients.

Ready to use?

If you already have your **doctor account**, then you can start.

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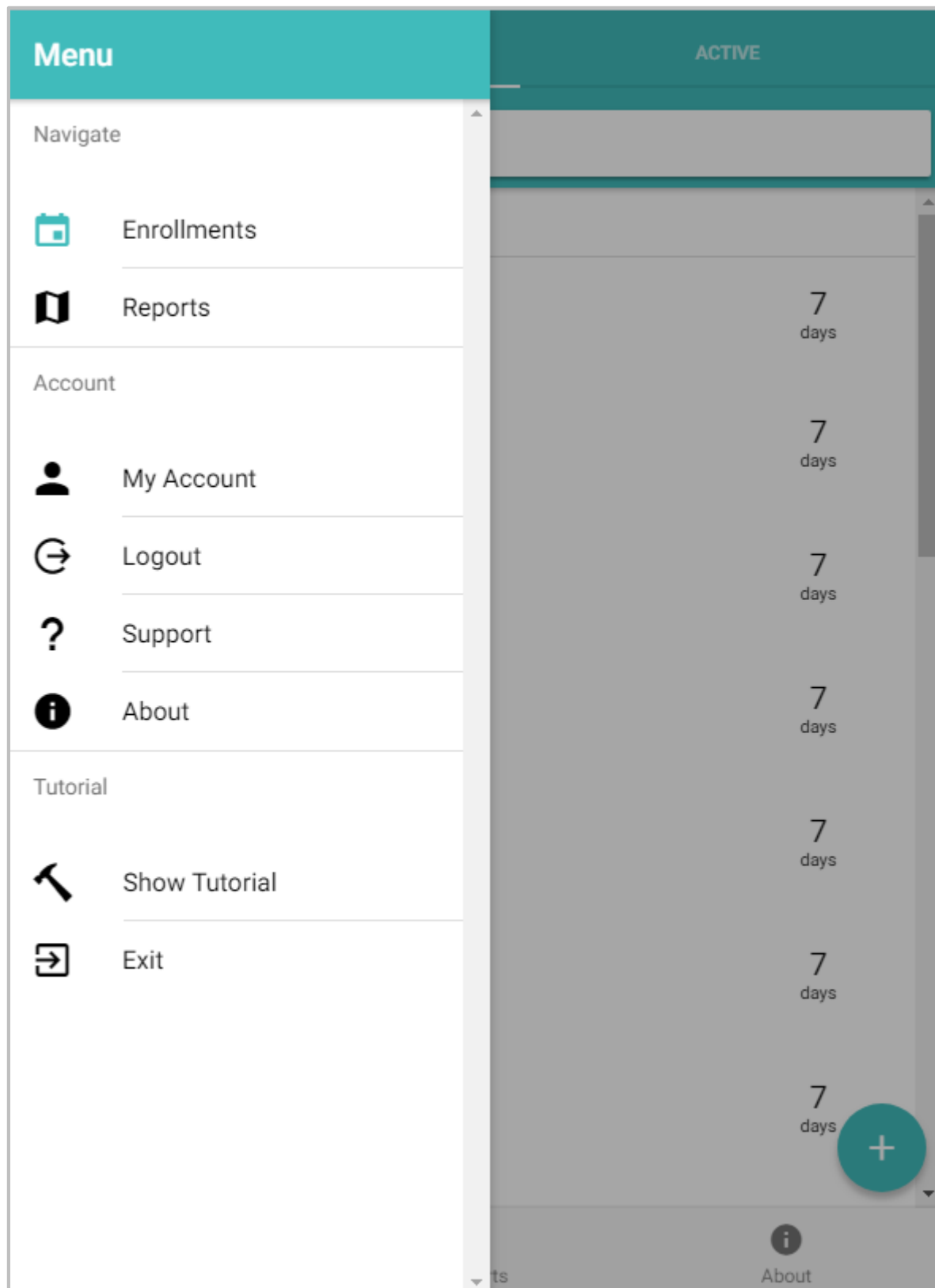
### Metadata

💬	Comments	null
🚴	Activities	null
📄	Symptoms	null

## ← 12309: OKTAY ZORLU

### Report Details

Status

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Report Date

29.01.2018 09:17

### Enrollment Details

Enrollment No

2428

Patient

OKTAY ZORLU

### Delivery Details

Delivery Status

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Manuel Gönderim .CallLogID: 83196ced-617d-45e1-8e2b-b52ec7e9a120, Konu: Mobil Kardiyak Telemetry İzlem Hizmet Sonu Raporu

### Metadata

Comments

null

Activities

null

Symptoms

null

SHOW REPORT

## My Account

tursoft

[Update Picture](#)

---

[Change Password](#)

---

[Support](#)

---

[Logout](#)

## Support

**ATTP**  
AVICENNA TELEMEDICINE PLATFORM

Enter your support message below

|

//

**SUBMIT**

About

**ATTP**  
AVICENNA TELEMEDICINE PLATFORM

## ATTP Mobile v0.9

ATTP Mobile is a full featured platform for Mobile Cardiac Management.

ERC Group, Engineering,  
Consultancy and R&D Services

[OPEN WEB PORTAL](#)

Enrollments

Reports

About